


stryker

SharpShooter


Inside-Out
Meniscal
Repair
System

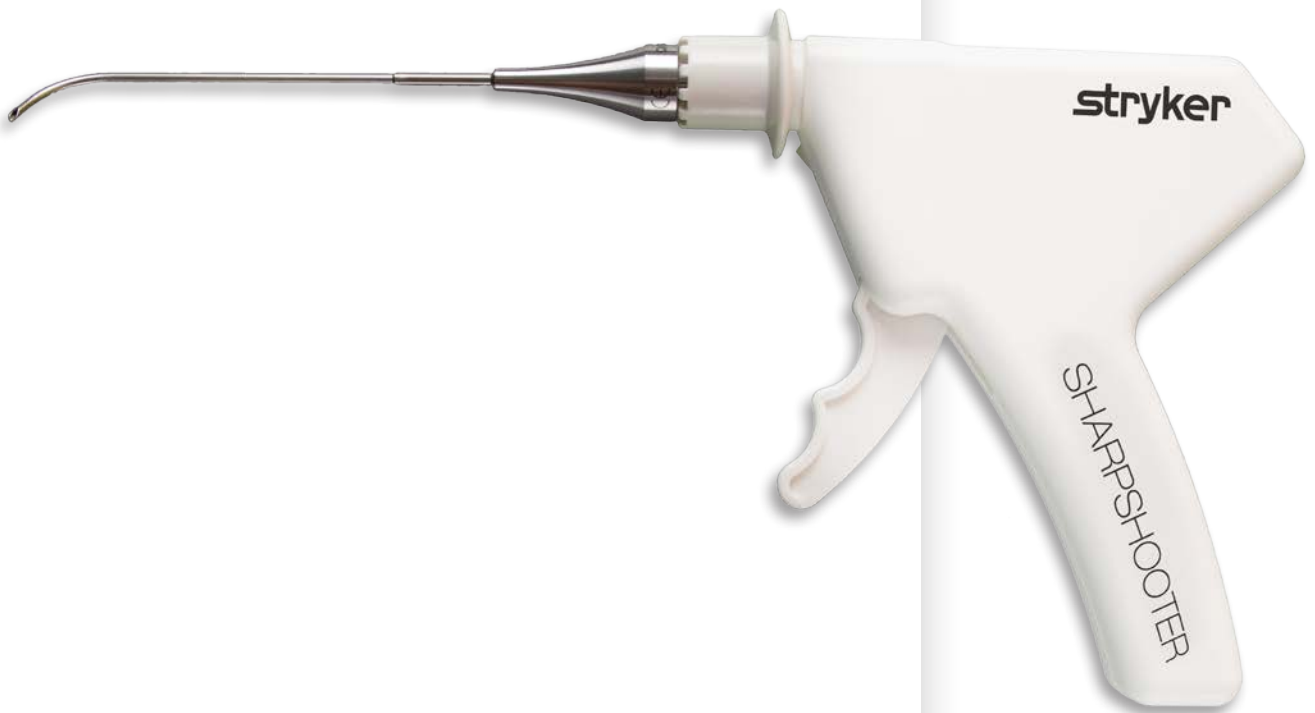

SharpShooter

Stryker's Inside-Out Meniscal Repair System

A simple, easy-to-use method for repairing the meniscus utilising an inside-out suture technique. The system includes a method to deliver either high strength 2-0 UHMWPE suture or 2-0 polyester braided suture using a patented delivery handle and six anatomic cannulas to reach all areas of the meniscus arthroscopically. The SharpShooter is designed to give the surgeon control of needle placement and delivery.

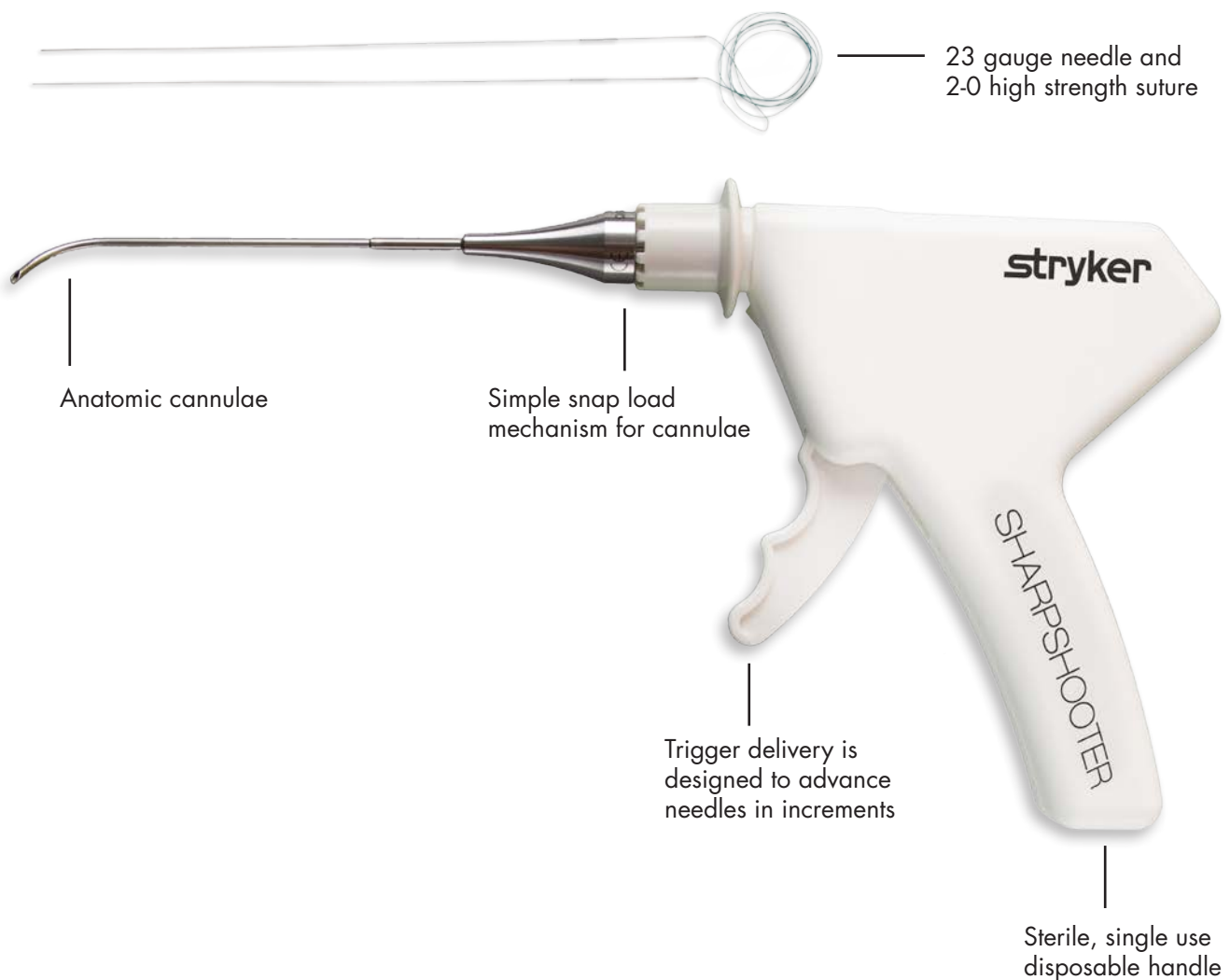

Features and benefits

Targeted Delivery

- SharpShooter Handle – A disposable suture passer that incrementally advances the SharpShooter needle
- SharpShooter Cannulas – Specifically designed cannulas that attach to the SharpShooter handle, designed to allow surgeons to gain access to all zones of the meniscus
 - Right or left options: posterior, middle, and anterior
- Designed to provide control of targeting needles

High-strength suture

- Size 2-0 UHMWPE suture or 2-0 polyester braided suture provides high suture strength demanded by surgeons

Efficient Design

- Single-handed needle advancement requires less operating assistance
- Double-armed, trocar-tipped suture needle incorporates visual markings designed to distinguish the distance traveled through the cannula before engaging the meniscus



Part Number	Description
4700	SharpShooter Handle
4701	SharpShooter Surgical Suture, braided PE
HS4704	SharpShooter Surgical Suture, UHMWPE
4750	SharpShooter Cannula Right Posterior
4751	SharpShooter Cannula Left Posterior
4752	SharpShooter Cannula Right Middle
4753	SharpShooter Cannula Left Middle
4754	SharpShooter Cannula Right Anterior
4755	SharpShooter Cannula Left Anterior
233050115	Popliteal Retractor Small
233050116	Popliteal Retractor Large

Sports Medicine

This document is intended solely for the use of healthcare professionals.

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product.

Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Sharpshooter, Stryker. All other trademarks are trademarks of their respective owners or holders.

The products depicted are CE marked in accordance with applicable EU Regulations and Directives.

This material is not intended for distribution outside the EU and EFTA.



1000902220EN
1000902220EN Rev B
DFS 04/2017
2017-13822

Stryker European Operations BV
Herikerbergweg 110
1101 CM Amsterdam
The Netherlands
t. +31 20 219 2600